



Protrieve™ Sheath


**Confidence During Complex DVT*
and IVC Thrombectomy Cases**



*Defined as bilateral DVT with or without thrombosed IVC


Trap Emboli

The **Protrieve Sheath** features a wall apposing funnel **designed to trap emboli** during complex DVT or IVC thrombectomy procedures

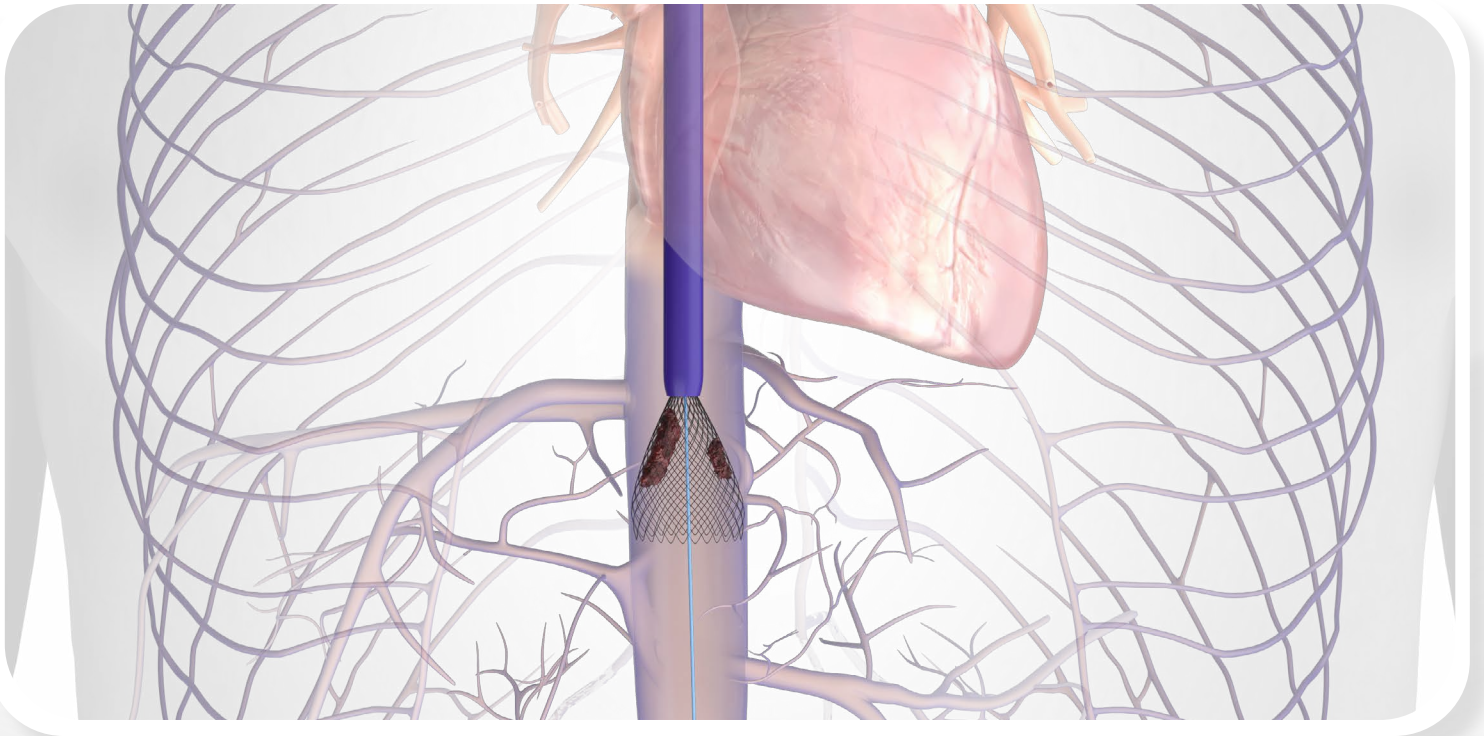


Peace of Mind

33.5mm self-expanding nitinol mesh funnel allows **full wall apposition in the IVC**

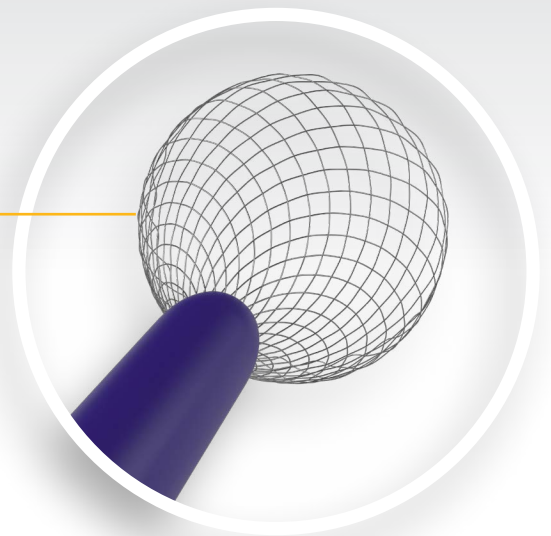

Seamless Venous Access


20F Sheath designed for **optimal access and positioning within the IVC**

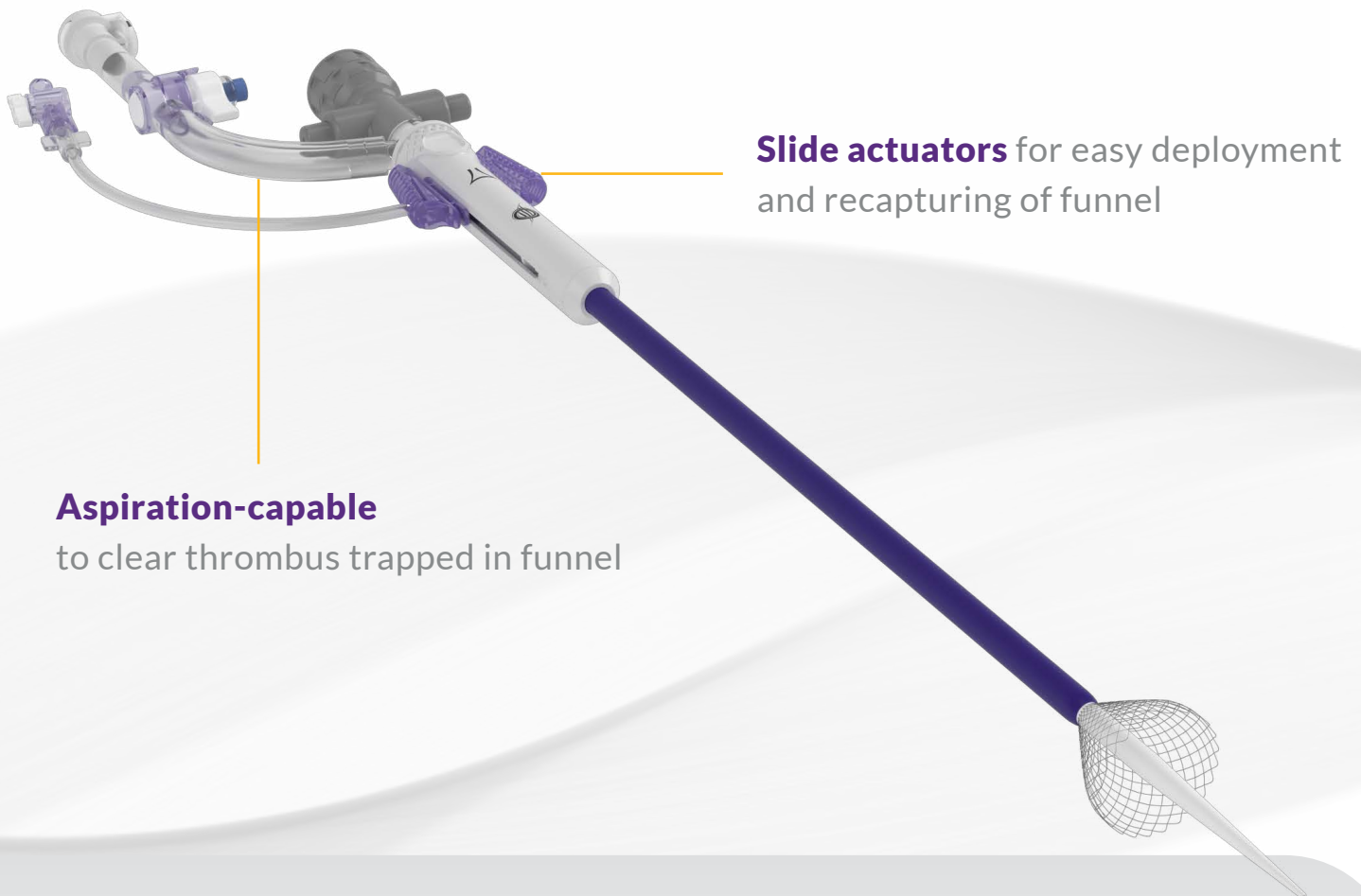


Slide actuators for easy deployment and recapturing of funnel

Aspiration-capable
to clear thrombus trapped in funnel

Compatible with **ClotTriever®** and **FlowTriever®** platforms[†]
for **simplification of complex thrombectomy**

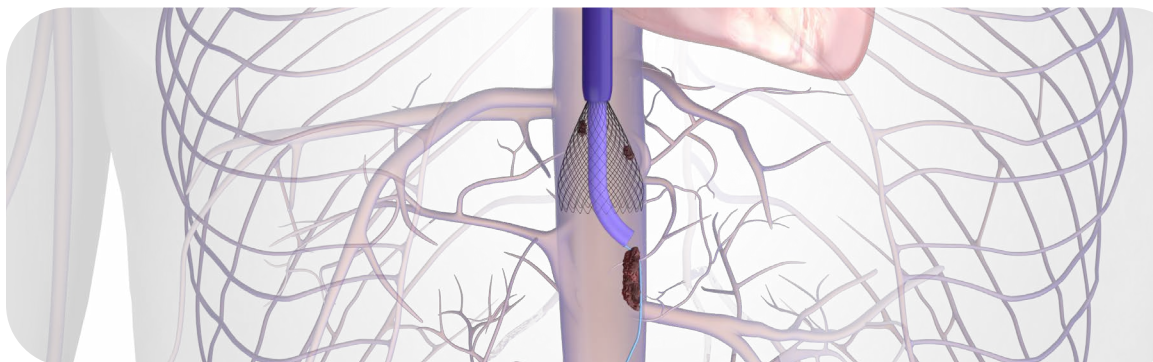

[†] Triever16®, Triever20® & Triever20 Curve®, FlowTriever catheters (S - M - L - XL), and FlowTriever2® catheter

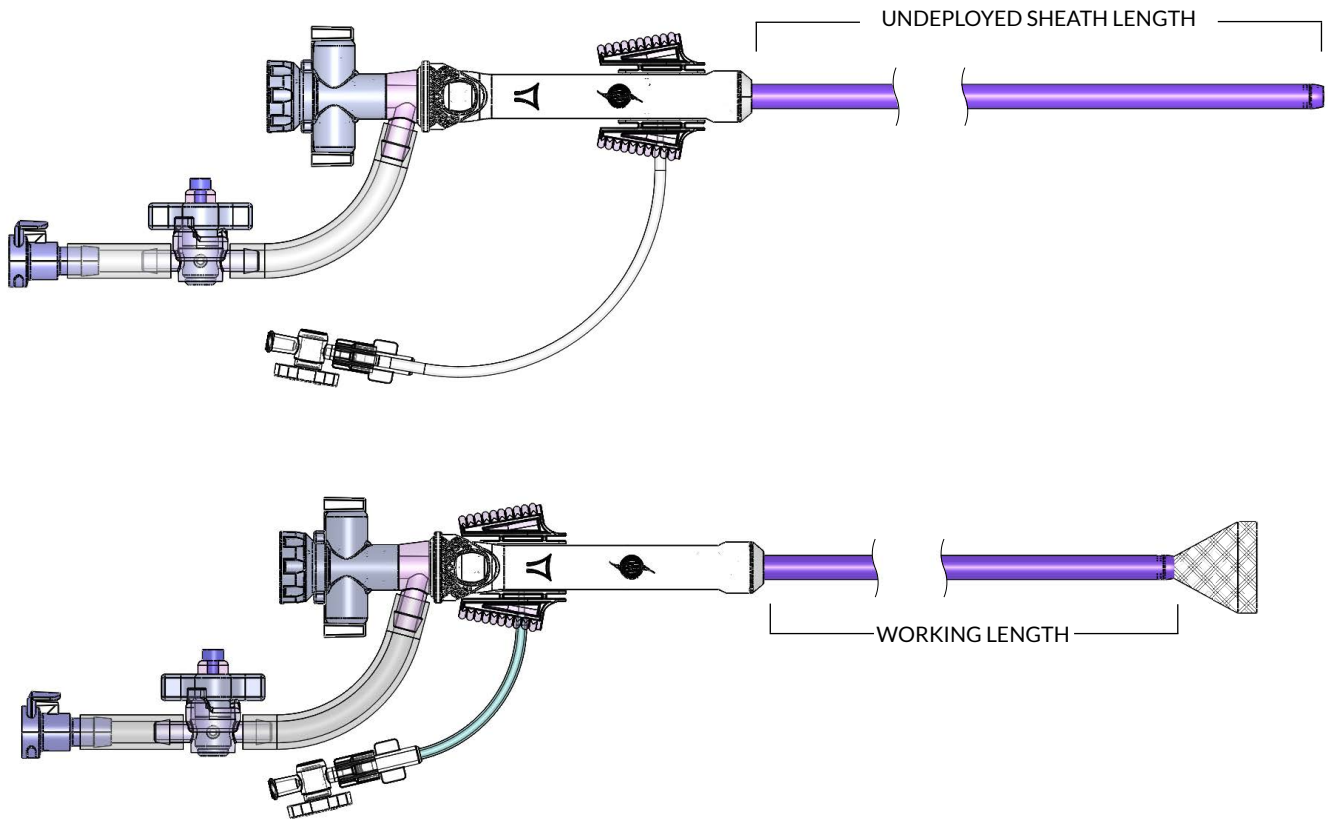

Model	Undeployed Sheath Length	Working Length	Outer Diameter	Inner Diameter	Deployed Funnel Diameter	Treatment Range
60-101	38 cm	32 cm	26F (8.8mm)	20F (6.8mm)	33.5 mm	17 -30 mm

Indications For Use: The ClotTriever thrombectomy system is indicated for: 1) The non-surgical removal of thrombi and emboli from blood vessels. 2) Injection, infusion, and/or aspiration of contrast media and other fluids into or from a blood vessel. The ClotTriever thrombectomy system is intended for use in the peripheral vasculature including deep vein thrombosis (DVT).

The FlowTriever retrieval/aspiration system is indicated for (1) the non-surgical removal of emboli and thrombi from blood vessels (2) injection, infusion and/or aspiration of contrast media and other fluids into or from a blood vessel. The FlowTriever retrieval/aspiration system is intended for use in the peripheral vasculature and for the treatment of pulmonary embolism. Triever catheters are intended for use in peripheral vasculature and for the treatment of pulmonary embolism. Triever catheters are also intended for use in treating clot in transit in the right atrium, but not in conjunction with FlowTriever Catheter. The FlowTriever2® catheter is indicated for the non-surgical removal of emboli and thrombi from peripheral blood vessels. Injection, infusion, and/or aspiration of contrast media and other fluids into or from a blood vessel. The FlowTriever2 catheter is intended for use in the peripheral vasculature.

See Instructions for Use for complete Indications for Use, contraindications, warnings, and precautions.

Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.

All trademarks are property of their respective owners.

MM-01314_Rev. A_EN_2022-09-21

